






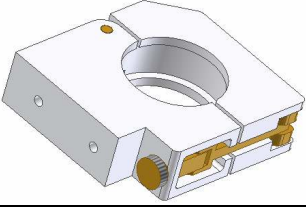
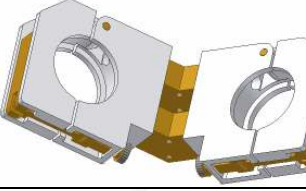
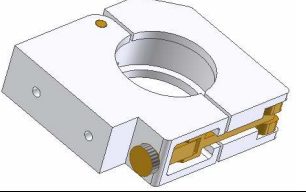







Accessories list for olfactometer and gustometer

Picture	name	Description	Fit for type	Order #
	Breathing sensor	Detects inhaling or exhaling of a subjects and initiate a TTL signal.	OL022/ OL023/ OL024/ OG001	LA-07-04005
	Feedback box	Allows to detect 2 feedback signals via 2 buttons. The software detects responds time from trigger release to keypress. The integrated lamp can be used for subject's preparation. It is connected to trigger opt. Works only in combination with Signal link box.	OL022/ OL023/ OL024/ OG001/ GU002	LA-07-02093
	Lamp	Allows to present a signal with trigger opt (for monorhinal use only!). Works only in combination with Signal link box.	OL022/ OL023/ OL024/ OG001/ GU002	LA-07-02094
	Signal link box	The triggerbox converts the 25-pin connection from all Burghart devices to a trigger-in, trigger-opt and trigger-out, as well as a connection for the Lamp, a feedback box and a USB Port for clear text transmission. At birhinal mode trigger out is trigger L(left) and trigger opt is trigger R (right)	OL022/ OL023/ OL024/ OG001/ GU002	LA-07-03290

	<p>Examination chair with holder</p>	<p>Complete examination chair with mounted holder for olfactometer connection hose and table for placing Burghart EEG system. Holder for hose must be order individually for the type.</p>	<p>OL022/ OL023/ OL024/ OG001</p>	<p>LA-07-02408</p>
	<p>Stand with holder</p>	<p>Stand with holder for olfactometer connection hose and table for placing Burghart EEG system which can be placed next to an existing chair. Holder for hose must be order individually for the type.</p>	<p>OL022/ OL023/ OL024/ OG001</p>	<p>LA-07-03359</p>
	<p>MRT holder</p>	<p>Holder which can be fixed at flatbed if fMRI scanner to hold hoses of olfactometer. Monorhinal and birhinal version available for OL023</p>	<p>OL023/OL022</p>	<p>LA-07-03349</p>

	holder for OL023 connection hose	Holder is intended to be used with Stand or with examination chair but can also be used for individual solutions.	OL023	LA-07-02505
	birhinal holder for OL023 connection hose	Holder is intended to be used with Stand or with examination chair but can also be used for individual solutions.	OL023	LA-07-00163
	holder for OL024 connection hose	Holder is intended to be used with Stand or with examination chair but can also be used for individual solutions.	OL024	LA-07-04036
	Humidification cartridge large	Module for fill amounts up to 50 ml. Originally developed for humidifying Dilution and Control flow. Can also be used in Odorant lines. Notice: Drop catchers are made of PE. This material might swell with some kind of odorants. Can be equipped with stainless steel drop catchers - please refer module spare parts list.	OL022/ OL023/ OL024	LA-07-04136
	Humidification cartridge small	Module for fill amounts up to 10 ml. Originally developed for odorant channels to save expensive odorant. It comes with stainless steel drop catchers.	OL022/ OL023/ OL024	LA-07-04135

	<p>cartridge for Pads</p>	<p>Module with 3 or more layers of filter paper which are clamped between meshes. Build for dropping a liquid on papers for odorization. Standard come with holder made of PE, but can be equipped with stainless steel drop catchers - please refer module spare parts list.</p>	<p>OL022/ OL023/ OL024</p>	<p>LA-07-04286</p>
	<p>cartridge for solids</p>	<p>Module with large space for solids like e.g. cotton pieces with sweat or potato chips. A filter paper between two meshes shall avoid large particles transported to the subjects nose. Can be equipped with stainless steel drop catchers - please refer module spare parts list.</p>	<p>OL022/ OL023/ OL024</p>	<p>LA-07-04116</p>
	<p>Burghart EEG 8channel</p>	<p>Especially build for recording event related evoked potentials. 8 channel device with complete recording and evaluation software. Presetting of filters for optimal results.</p>	<p>OL022/ OL023/ OL024/ OG001/ GU002</p>	<p>LA-07-00124</p>

all pictures and descriptions are subject to technical alternations and might differ from delivered devices